




Rondish HelpSignal Emergency Dialer

NEW!

Remote supervision of loved ones, without monthly charges

- If you or your professional caregiver are looking after a parent, friend or someone you are responsible for, and can't be with them all the time, HelpSignal is here to help you.
- The plug and play HelpSignal uses any home telephone socket and provides the ability to remotely monitor someone, enabling them to call you quickly when needed.
- HelpSignal consists of an autodialer which is triggered using a wireless pendant worn by the user, at the simple press of a button.
- It automatically dials up to 4 numbers of relatives, friends, neighbors or professional caregivers who are remotely supervising the person under care.
- Once connected, a pre-recorded message giving caller information is played, after which it switches to speakerphone function allowing the carer to talk to the user and assess the nature of the caller's condition.
- HelpSignal may also be triggered from the main unit by pressing any one of the 3 orange "Emergency Call" buttons below the dial pad.
 - The self-recorded short help message, repeats 3 times (maximum 20 seconds).
 - Each of 4 designated emergency numbers will be dialed up to 15 times maximum, or a minimum of 2 successful calls to each number, unless "Reset" is pressed by the user. This ensures the greatest chance of calls being received and responded to by carers.

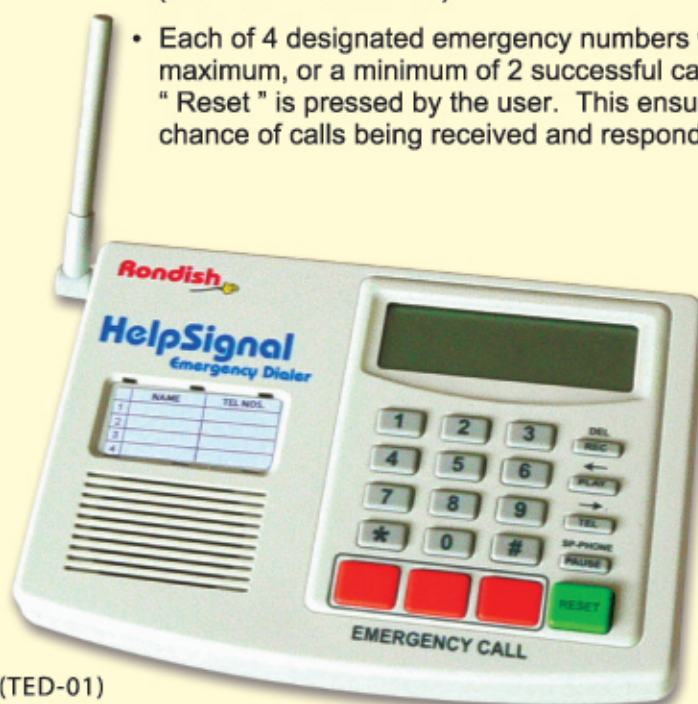

(TED-01)



(TTX-01)


HelpSignal Emergency Dialer

Standard Features

- Attractive, small, HelpSignal main unit with a clear LCD for ease of use.
- The HelpSignal supports up to 7 uniquely-coded Wireless help buttons, so 6 additional users may be easily logged onto the same main unit, dialing the same number. Flashing LED on wireless button confirms activation.
- Designed to be easy to program and set up - plug and play!
- Comes with battery back-up, resident's phone override in an emergency and reassurance tones during dialing so the user knows progress is being made to get help.
- The HelpSignal Emergency Dialer may also be used as a normal speakerphone. Just press the SP-PHONE button to connect the line and dial a telephone number on the keypad.

Specifications

| | |
|-----------------------------------|---|
| Time interval between each redial | 40 seconds between each call until answer |
| Dialing | DTMF Tone dialing |
| Maximum times of redial | 15 attempts |
| Duration of recorded message | Maximum 20 seconds |
| Wireless range (HelpSignal) | 50M (315 MHz - USA, 433 MHz - Europe) |
| Battery use | 4 x AA (Dialer), 1 x 23A (Wireless help button) |
| AC adapter | 9V 300mA |
| External siren output | 5 second pulse, 100mA to ground |
| External trigger input | Volt free open to close contact |

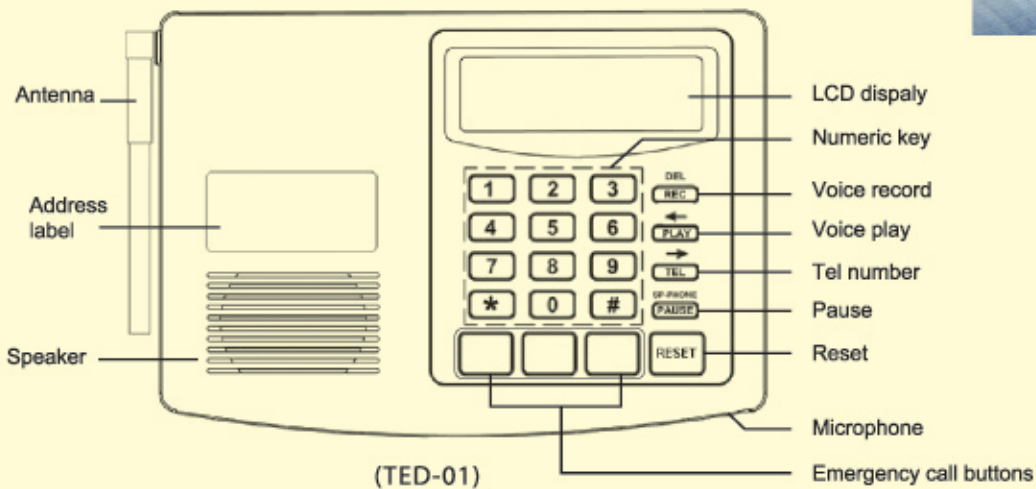

NOTE: HelpSignal has different country versions, primarily European and US types.

Distributed By:



Hong Kong Office:
Rondish Company Limited
Tel: +852 2545 2057
Fax: +852 2541 7411
info@rondish.com
www.rondish.com



USA Office:
RondishUSA LLC
Tel: +1 858 729 0251
Fax: +1 815 346 5396
sales@rondishUSA.com
www.rondishUSA.com